Blood Agar Base

REF CM0055B and CM0055T

EN

Intended Use

Blood Agar Base (CM0055B and CM0055T) devices are a nonselective, general-purpose medium intended to be used with the addition of blood for the isolation and culture of pathogenic and non-pathogenic microorganisms from clinical samples (e.g. cerebrospinal fluid (CSF), pus and wound swabs, sputum, eye swabs and ear, nose and throat swabs).

Blood Agar Base (CM0055B and CM0055T) devices are used in a diagnostic workflow to aid clinicians in determining potential treatment options for patients suspected of having a microbial infection.

The devices are for professional use only, are not automated and nor are they companion diagnostics.

Summary and Explanation

Non-fastidious microorganisms are organisms that grow and replicate without special nutritional supplements or conditions. They can synthesize all organic molecules required for their growth. Non-fastidious microorganisms can grow in nutrient sparse environments and can survive in a wider temperature range. This group of bacteria are involved in causing illnesses and infections in humans and includes *Escherichia coli* and *Staphylococcus aureus*^{1,2,3,4,5}.

Fastidious organisms are organisms that have complex or particular nutritional requirements in order to grow. Fastidious bacteria are not always easy to grow in culture and may not be easily obtained from a sample. This can be problematic in the laboratory, where a specific culture medium is required for growth.

Principle of Method

Blood Agar Base (CM0055B and CM0055T) devices are a nonselective, general-purpose medium which may be enriched with blood to facilitate the growth of a wide variety of microorganisms, including fastidious organisms, from clinical specimens.

Without any additions, the medium can be employed as a nutrient agar used, for example, for short-term maintenance of stock cultures. With added blood, the medium is suitable for cultivation of many fastidious organisms. The addition of blood allows the medium to become suitable for the determination of haemolytic reactions which is an important diagnostic criterion for many organisms.

To encourage growth of microorganisms tryptone and peptones, which provide nitrogen, carbon and amino acids, and yeast extract, which provides growth promoting B vitamins, are included. Sodium chloride maintains the osmotic balance. Agar is the solidifying agent.

Typical Formula

	<u>grams per litre</u>
"Lab-Lemco" powder	10.0
Peptone	10.0
Sodium chloride	5.0
Agar	15.0

Materials Provided

CM0055B: 500g of Blood Agar Base powder that yields approximately 12.5L after reconstitution.

CM0055T: 5kg of Blood Agar Base powder that yields approximately 125L after reconstitution.

Materials Required but Not Supplied

- Inoculating loops, swabs, collection containers
- Incubators
- Quality control organisms
- Optional Selective supplements (SR0050B/C/D/E/F/R)
- Petri dish

Storage

- Store product in its original packaging between 10°C and 30°C.
- Keep container tightly closed.
- The product may be used until the expiry date stated on the label.
- Protect from moisture.
- Store away from light.
- Allow reconstituted product to equilibrate to room temperature before use.

Once reconstituted, store media between 2°C and 10°C.

Warnings and Precautions

- Do not inhale. May cause allergy or asthma symptoms or difficulty breathing if inhaled.
- Causes serious eye irritation.
- May cause an allergic skin reaction.
- If on skin wash with plenty of soap and water.
- If in eyes, rinse cautiously with water for several minutes.
- Remove contact lenses, if present and easy to do. Continue rinsing. If eye irritation persists, seek medical advice/attention.
- If inhaled, if breathing is difficult, remove subject to fresh air and keep in a position comfortable for breathing. If experiencing respiratory symptoms, call a POISON CENTER or doctor/physician.
- For in vitro diagnostic use only.
- For professional use only.
- Inspect the product packaging before first use.
- Do not use the product if there is any visible damage to the packaging (pot or cap).
- Do not use the product beyond the stated expiry date.
- Do not use the device if signs of contamination are present.
- It is the responsibility of each laboratory to manage waste produced according to their nature and degree of hazard and to have them treated or disposed of in accordance with any federal, state and local applicable regulations. Directions should be read and followed carefully. This includes the disposal of used or unused reagents as well as any other contaminated disposable material following procedures for infectious or potentially infectious products.
- Ensure the lid of the container is kept tightly closed after first opening and between use to minimise moisture ingress, which may result in incorrect product performance.

Refer to the Safety Data Sheet (SDS) for safe handling and disposal of the product (www.thermofisher.com).

Serious Incidents

Any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the relevant regulatory authority in which the user and/or the patient is established.

Specimen Collection, Handling and Storage

Specimen should be collected and handled following local recommended guidelines, such as the UK Standards for Microbiology Investigations (UK SMI) B 29 (Public Health England, 2020).

Procedure

Suspend 40g in 1 litre of distilled water. Bring to the boil to dissolve completely. Sterilize by autoclaving at 121°C for 15 minutes. Cool to 50°C. Mix well and pour into sterile Petri dishes. For blood agar, enrich with 7% v/v sterile defibrinated blood (SR0050B/C/D/E/F/R).

Quality Control

It is the responsibility of the user to perform Quality Control testing taking into account the intended use of the medium, and in accordance with any local applicable regulations (frequency, number of strains, incubation temperature etc.).

The performance of this medium can be verified by testing the following reference strains.

Incubation Conditions: 24 h @ 37°C aerobic

Positive Controls	
Colony count is $\geq 70\%$ of the control medium count	
<i>Staphylococcus aureus</i> ATCC® 25923	1-1.5mm straw/white colonies
<i>Pseudomonas aeruginosa</i> ATCC® 27853	1-2mm straw colonies with or without green pigmentation
<i>Escherichia coli</i> ATCC®25922	1-2mm straw colonies
Positive Controls enriched with 7% v/v horse blood	
Colony count is $\geq 70\%$ of the control medium count	
<i>Streptococcus pyogenes</i> ATCC® 19615	Pinpoint-0.75mm pale colonies, β haemolysis
<i>Streptococcus pneumoniae</i> ATCC® 6303	1-1.5mm grey/green mucoid colonies, α haemolysis
<i>Streptococcus pneumoniae</i> ATCC® 6305	0.25-0.5mm grey/green colonies, α haemolysis
<i>Haemophilus influenzae</i> ATCC® 19418	Pinpoint-0.25mm pale colonies

Zones of growth/no growth surrounding X, V and X+V factor discs when plain plates are inoculated with the following organisms and incubated at 37°C for 18 hours:			
	X	V	X+V
<i>Haemophilus influenzae</i> ATCC® 49247™	0	0	≥15 mm
<i>Haemophilus parainfluenzae</i> ATCC® 33392™	0	≥20 mm	≥20 mm

Zones of inhibition with Bacitracin discs when inoculated on 7% v/v horse blood plates Incubation conditions: 37°C for 18 hours.	
	Zones of inhibition
<i>Streptococcus pyogenes</i> ATCC® 19615™	10-20mm

Limitations

Blood Agar Base (CM0055B and CM0055T) devices are a nonselective, general-purpose medium intended to be used with the addition of blood for the isolation and culture of pathogenic and non-pathogenic microorganisms from clinical samples (e.g. cerebrospinal fluid (CSF), pus and wound swabs, sputum, eye swabs and ear, nose and throat swabs). As with all media, atypical organisms may give anomalous reactions. A small number of atypical strains may give a weak reaction and fail to grow, especially when low numbers are present in the sample. Further biochemical and serological tests must be carried out for further identification.

Performance Characteristics

Accuracy has been demonstrated through review of the QC data. Correct isolation and culture of pathogenic and non-pathogenic microorganisms from clinical samples is confirmed by the inclusion of a well-characterised isolate in the QC processes performed as part of the manufacture of each batch of the devices. The precision of Blood Agar Base (CM0055B and CM0055T) was demonstrated by an overall pass rate of 100% obtained for the product over 10 months of testing (30.JUL.2021 – 21.MAY.2022 10 batches). This shows that the performance is reproducible.

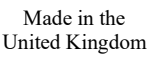
These devices have been tested in-house as part of the QC process since they were first manufactured in 1959. When using 10 - 100 cfu inoculum of *Escherichia coli*, *Pseudomonas aeruginosa*, *Staphylococcus aureus*, *Streptococcus pyogenes*, *Streptococcus pneumoniae* and *Haemophilus influenzae* and incubating the device at 37°C for 24 hours, the user can recover organisms with colony size and morphology as listed in this document

Bibliography

1. World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Accessed 29.03.2022
4. Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

Symbol Legend

Symbol	Definition
	Catalogue number
	In Vitro Diagnostic Medical Device
	Batch code
	Temperature limit

	Use-by date
	Keep away from sunlight
	Do not re-use
	Consult instructions for use or consult electronic instructions for use
	Contains sufficient for <n> tests
	Do not use if packaging damaged and consult instructions for use
	Manufacturer
	Authorized representative in the European Community/ European Union
	European Conformity Assessment
	UK Conformity Assessment
	Unique device identifier
	Importer - To indicate the entity importing the medical device into the locale. Applicable to the European Union
	Made in the United Kingdom

ATCC Licensed
Derivative

©2022 Thermo Fisher Scientific Inc. All rights reserved. ATCC and ATCC catalogue marks are a trademark of American Type Culture Collection.
All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

Oxoid Limited, Wade Road, Basingstoke, Hampshire, RG24 8PW, UK

For technical assistance please contact your local distributor.

Revision information

Version	Date of issue and modifications introduced
2.0	2023-12-12

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Litteratur

1. Verdenssundhedsorganisationen (WHO). 2017. "Microbial Fact Sheets: Guidelines for drinking-water quality". https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. "Identification of Enterobacteriaceae". UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. "UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens." Tilgængelig på: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Tilgængeligt 29.03.2022
4. Public Health England. 2020a. "Identification of Staphylococcus species, Micrococcus species and Rothia species." UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. "Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners." Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Literaturverzeichnis

1. Weltgesundheitsorganisation (WHO). 2017. „Microbial Fact Sheets: Guidelines for drinking-water quality“. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. „Identification of Enterobacteriaceae“. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. „UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.“ Verfügbar unter: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Zugriff am 29.03.2022
4. Public Health England. 2020a. „Identification of Staphylococcus species, Micrococcus species and Rothia species“. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. „Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.“ Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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1. World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Accessed 29.03.2022
4. Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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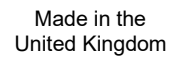
- Organización Mundial de la Salud (OMS). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
- Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
- Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Consultado el 29/03/2022
- Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
- Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Oxoid Limited, Wade Road, Basingstoke, Hampshire, RG24 8PW, Reino Unido

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Kirjandus

1. Maailma Terviseorganisatsioon (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Külastatud 29.03.2022
4. Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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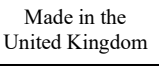
Bibliographie

1. Organisation mondiale de la santé (OMS). 2017. «Microbial Fact Sheets: Guidelines for drinking-water quality». https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. «Identification of Enterobacteriaceae». UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015 b. «UK SMI B57: Investigation of bronchoalveolar lavage, expectorations and Associateschants». Disponible à l'adresse : <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Consulté le 29/03/2022
4. Public Health England. 2020a. «Identification of Staphylococcus species, Micrococcus species and Rothia species». UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. «Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners». Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Oxoid Limited, Wade Road, Basingstoke, Hampshire, RG24 8PW, Royaume-Uni

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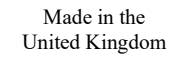
Bibliografia

1. World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. "Identification of Enterobacteriaceae". UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Accessed 29.03.2022
4. Public Health England. 2020a. "Identification of Staphylococcus species, Micrococcus species and Rotaia species". UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. "Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners". Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Oxoid Limited, Wade Road, Basingstoke, Hampshire, RG24 8PW, Regno Unito

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- World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
- Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
- Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Pieejams vietnē: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Piekļuve veikta 29.03.2022.
- Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
- Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Bibliografi

1. World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. 'Identification of Enterobacteriaceae.' UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Accessed 29.03.2022.
4. Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species.' UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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- Światowa Organizacja Zdrowia (WHO). 2017. „Microbial Fact Sheets: Guidelines for drinking-water quality”. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
- Public Health England. 2015a. „Identification of Enterobacteriaceae”. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
- Public Health England. 2015b. „UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.” [online:] <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Dostęp: 29.03.2022 r.
- Public Health England. 2020a. „Identification of Staphylococcus species, Micrococcus species and Rothia species”. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
- Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd i wsp., 2006. „Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.” Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C–24C. <https://doi.org/10.1155/2006/971352>.

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Bibliografia

1. World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Acedido em 29-Mar-2022
4. Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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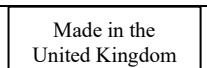
Bibliografie

1. Organizația Mondială a Sănătății (OMS). 2017. „Microbial Fact Sheets: Guidelines for drinking-water quality”. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. „Identification of Enterobacteriaceae”. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. „UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.” Disponibilă la: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Accesat pe 29.03.2022.
4. Public Health England. 2020a. „Identification of Staphylococcus species, Micrococcus species and Rothia species”. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd, et al. 2006. „Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.” Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Oxoid Limited, Wade Road, Basingstoke, Hampshire, RG24 8PW, Regatul Unit

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Bibliografia

1. World Health Organization (WHO). 2017. 'Microbial Fact Sheets: Guidelines for drinking-water quality'. https://www.who.int/water_sanitation_health/publications/gdwq4-with-add1-chap11.pdf.
2. Public Health England. 2015a. 'Identification of Enterobacteriaceae'. UK Standards for Microbiology Investigations ID 16 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/423601/ID_16i4.pdf.
3. Public Health England. 2015b. 'UK SMI B57: Investigation of bronchoalveolar lavage, sputum and associated specimens.' Available at: <https://www.gov.uk/government/publications/smi-b-57-investigation-of-bronchoalveolar-lavage-sputum-and-associated-specimens>. Accessed 29.03.2022
4. Public Health England. 2020a. 'Identification of Staphylococcus species, Micrococcus species and Rothia species'. UK Standards for Microbiology Investigations ID 07 (4). https://assets.publishing.service.gov.uk/government/uploads/system/uploads/attachment_data/file/887570/UK_SMI_ID_07i4.pdf.
5. Barton, Michelle, Michael Hawkes, Dorothy Moore, John Conly, Lindsay Nicolle, Upton Allen, Nora Boyd a kol. 2006. 'Guidelines for the Prevention and Management of Community-Associated Methicillin-Resistant Staphylococcus aureus: A Perspective for Canadian Health Care Practitioners.' Canadian Journal of Infectious Diseases and Medical Microbiology 17 (c): 4C-24C. <https://doi.org/10.1155/2006/971352>.

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Všetky ostatné ochranné známky sú vlastníctvom Thermo Fisher Scientific Inc. a jej dcérskych spoločností.
Oxoid Limited, Wade Road, Basingstoke, Hampshire, RG24 8PW, SPOJENÉ KRÁĽOVSTVO