Columbia Blood Agar Base

EN

REF CM0331B, CM0331K, CM0331M, CM0331R, CM0331T, CM0331U CM0331W and CM0331Y

Intended Use

Columbia Blood Agar Base devices are nonselective, general-purpose medium which may be enriched with blood or supplements intended for the isolation and culture of fastidious microorganisms from clinical samples (e.g., wounds, throat, genital, nose, skin, etc.).

Columbia Blood Agar Base devices are used in a diagnostic workflow to aid clinicians in determining potential treatment options for patients suspected of having bacterial infections.

These devices are for professional use only, is not automated, nor is it a companion diagnostic.

Summary and Explanation

Gram-positive and Gram-negative microorganisms exist in a wide range of ecological niches, including environmental bacteria, commensals, and human pathogens¹. Prominent Gram-positive human pathogens include *Streptococcus* species, *Staphylococcus* species². The *Streptococcus* genus is comprised of many pathogenic commensals and some of the important *Streptococcus* species are *Streptococcus pneumoniae*, *Streptococcus viridans* and *Streptococcus pyogenes*. and the optimum growth temperature of these species is identified as 37°C¹.

Other prominent Gram-positive human pathogens include *Staphylococcus aureus*³. *S. aureus* is an opportunistic pathogen, entering the body through cuts or abrasions where it can cause a spectrum of disease in both immunocompetent and immunocompromised hosts and can range from superficial skin infections such as cellulitis and abscesses to life-threatening illnesses such as pneumonia, toxic shock syndrome and sepsis. *S. aureus* is one of the leading causes of hospital acquired infections (HAIs), further complicated by the rapid emergence of antibiotic resistance in this pathogen³. Gram-negative pathogen that are involved in blood stream related infections (BSI) in patients receiving invasive operations are *Escherichia coli* and is also considered as the second most common bacterium that is isolated from the blood of neonates identified with sepsis¹. *Neisseria* species is another human pathogen that is nutritionally fastidious and haemolytic in nature and can cause meningococcal disease in humans⁴. Among anaerobic pathogenic bacteria *Clostridium* belongs to Clostridiaceae family and a few strains are known to cause illness in humans⁵.

Principle of Method

Columbia Blood Agar Base devices contains peptones in the media, which provide essential growth factors such as nitrogen, carbon, vitamins and trace elements necessary for bacterial growth. Starch is added to absorb any toxic metabolites. Columbia Blood Agar Base devices when supplemented with defibrinated blood (horse/sheep) which acts as a nutrient source and allows the determination of the haemolytic properties of the isolate. The presence and pattern of haemolysis, the size of the zone of haemolysis and colony morphology can be used to aid in the identification of isolates.

Typical Formula

	grams per litre
Special peptone	23
Soluble starch	1
Sodium chloride	5
Agar	10

Materials Provided

CM0331B: 500g of Columbia Blood Agar powder that yields approximately 12.8L after reconstitution.

CM0331K: 25kg of Columbia Blood Agar powder that yields approximately 641L after reconstitution.

CM0331M: 3 x 2.34kg of Columbia Blood Agar powder

CM0331R: 2.5kg of Columbia Blood Agar powder that yields approximately 64.1L after reconstitution.

CM0331T: 5kg of Columbia Blood Agar powder that yields approximately 128L after reconstitution.

CM0331U: 6 X 4.19g of Columbia Blood Agar powder

CM0331W: 3 X 3.9kg of Columbia Blood Agar powder

CM0331Y: 3 X 4.19kg of Columbia Blood Agar powder

Materials Required but Not Supplied

- Inoculating loops, swabs, collection containers
- Incubators
- Quality control organisms
- Petri dish

Optional Selective supplements

- Laked Horse Blood (SR0048C)
- Defibrinated Horse Blood (SR0050B/C/D/E/F/R)
- Defibrinated Sheep Blood (SR0051B/C/D/E/R)

- Campylobacter Selective Supplement (Skirrow) (SR0069E/H)
- Staph/Strep Selective Supplement (SR0070E)
- Campylobacter Selective Supplement(Butzler) (SR0085E)
- Campylobacter Selective Supplement (Blaser-Wang) (SR0098E)
- Streptococcus Selective Supplement (SR0126E)
- Helicobacter pylori Selective Supplement (DENT)(SR0147E)
- Modified CNA Selective Supplement (SR0176G)
- Campylobacter Growth Supplement (SR0232E)

Storage

- Store product in its original packaging between 10°C and 30°C.
- Keep container tightly closed.
- The product may be used until the expiry date stated on the label.
- Protect from moisture.
- Store away from light.
- Allow reconstituted product to equilibrate to room temperature before use.

Once reconstituted, store media between 2°C and 8°C.

Warnings and Precautions

- Do not inhale. May cause allergy or asthma symptoms or difficulty breathing if inhaled.
- Causes serious eye irritation.
- May cause an allergic skin reaction.
- If on skin wash with plenty of soap and water.
- If in eyes, rinse cautiously with water for several minutes.
- Remove contact lenses, if present and easy to do. Continue rinsing. If eye irritation persists, seek medical advice/attention.
- If inhaled, if breathing is difficult, remove subject to fresh air and keep in a position comfortable for breathing. If experiencing respiratory symptoms, call a POISON CENTER or doctor/physician.
- For in vitro diagnostic use only.
- For professional use only.
- Inspect the product packaging before first use.
- Do not use the product if there is any visible damage to the packaging (pot or cap).
- Do not use the product beyond the stated expiry date.
- Do not use the device if signs of contamination are present.
- It is the responsibility of each laboratory to manage waste produced according to their nature and degree of hazard and to have them treated or disposed of in accordance with any federal, state and local applicable regulations. Directions should be read and followed carefully. This includes the disposal of used or unused reagents as well as any other contaminated disposable material following procedures for infectious or potentially infectious products.
- Ensure the lid of the container is kept tightly closed after first opening and between use to minimise moisture ingress, which may result in incorrect product performance.

Refer to the Safety Data Sheet (SDS) for safe handling and disposal of the product (www.thermofisher.com).

Serious Incidents

Any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the relevant regulatory authority in which the user and/or the patient is established.

Specimen Collection, Handling and Storage

Specimen should be collected and handled following local recommended guidelines, such as the UK Standards for Microbiology Investigations (UK SMI) ID 4, ID 6, ID 7, ID 8, ID 12, ID 23 and Q5.

Procedure

Suspend 39g in 1 litre of distilled water. Bring to the boil to dissolve completely. Sterilize by autoclaving at 121°C for 15 minutes. Cool to 50°C. Mix well and pour into sterile Petri dishes.

For blood agar, enrich with 5% v/v sterile defibrinated blood.

Quality Control

It is the responsibility of the user to perform Quality Control testing taking into account the intended use of the medium, and in accordance with any local applicable regulations (frequency, number of strains, incubation temperature etc.).

The performance of this medium can be verified by testing the following reference strains.

Positive Controls Incubation Conditions: 24 h @ 37°C aerobic Inoculum level: 10-100 cfu Colony count is ≥ 70% of the control medium count	
<i>Staphylococcus aureus</i> ATCC® 25923™	1-2 mm cream colonies
Positive Controls Incubation Conditions: 48 h @ 37°C anaerobic Colony count is ≥ 70% of the control medium count	

<i>Clostridium sporogenes</i> ATCC® 19404™	1-2 mm pale straw colonies
<i>Clostridium sporogenes</i> ATCC® 11437™	1-2 mm pale straw colonies
Positive Controls Incubation Conditions: 48 h @ 37°C 5% v/v sheep blood plates microaerophilic Colony count is ≥ 70% of the control medium count	
<i>Neisseria gonorrhoeae</i> NCTC 11148	1-2mm grey/brown colonies
Positive Controls Incubation Conditions: 24 h @ 37°C aerobic and enriched with 5% v/v horse blood (SR0050) Colony count is ≥ 70% of the control medium count	
<i>Streptococcus pyogenes</i> ATCC® 19615™	0.25-1mm pale straw colonies, β haemolysis
<i>Streptococcus pneumoniae</i> ATCC® 6305™	0.5-1mm grey/green colonies, α haemolysis

Zones of growth/no growth surrounding X, V and X+V factor discs when plain plates are inoculated with the following organisms and incubated at 37°C for 18 hours:			
	X	V	X+V
<i>Haemophilus influenzae</i> ATCC® 49247™	0	0	≥15 mm
<i>Haemophilus parainfluenzae</i> ATCC® 33392™	0	≥20 mm	≥20 mm

Zones of inhibition with Bacitracin discs when inoculated on 7% v/v horse blood plates Incubation conditions: 37°C for 18 hours.	
	Zones of inhibition
<i>Streptococcus pyogenes</i> ATCC® 19615™	10-20mm

Testing performed in accordance with ISO11133:2014 Incubation conditions: 41.5 ± 2°C for 44 ± 4 hours, microaerophilic Inoculum level: 50-120 cfu Colony count is ≥ 70% of the control medium count	
<i>Campylobacter jejuni</i> ATCC® 29428™	1-3mm grey, mucoid colonies
<i>Campylobacter jejuni</i> ATCC® 33291™	1-3mm grey, mucoid colonies
<i>Campylobacter coli</i> ATCC® 43478™	1-3mm grey, mucoid colonies

Testing performed in accordance with current CLSI M22 A Incubation conditions: 35°C for 18-24 hours Inoculum level: 10 ² – 10 ⁴ cfu	
<i>Streptococcus pyogenes</i> ATCC® 19615™	(++ to +++) 0.25-1 mm pale straw colonies, with β (beta) haemolysis
<i>Streptococcus pneumoniae</i> ATCC® 6305™	(++ to +++) 0.5-1mm grey/green colonies, with α (alpha) haemolysis
<i>Staphylococcus aureus</i> ATCC® 25923™	(++ to +++) 1-2 mm cream colonies
<i>Escherichia coli</i> ATCC® 25922™	(++ to +++) 1-3 mm cream colonies

Limitations

As with all media, atypical organisms may give anomalous reactions. A small number of atypical strains may give a weak reaction

and fail to grow, especially when low numbers are present in the sample. Further biochemical and serological tests must be carried out for further identification.

Performance Characteristics

Accuracy has been demonstrated through review of the QC data. Correct detection of microorganism is confirmed by the inclusion of a well-characterised isolates in the QC processes performed as part of the manufacture of each batch of the device. The precision of Columbia Blood Agar Base (CM0331B) was demonstrated by an overall pass rate of 100% obtained for the product over a 10-day period (04.07.2022 – 14.07.2022; 10 batches). The precision of Columbia Blood Agar Base (CM0331K) was demonstrated by an overall pass rate of 100% in two months of testing (27.02.2012 – 20.04.2012; 2 batches). The precision of Columbia Blood Agar Base (CM0331R) was demonstrated by an overall pass rate of 100% in 10 years of testing (23.07.2012 – 04.07.2022; 10 batches). The precision of Columbia Blood Agar Base (CM0331T) was demonstrated by an overall pass rate of 100% in seven years of testing (13.05.2015 – 04.07.2022; 10 batches). The precision of Columbia Blood Agar Base (CM0331U) was demonstrated with an overall pass rate of 100% in eight months of testing (19.06.2017 – 15.02.2018; 10 batches). Columbia Blood Agar Base (CM0331W) has demonstrated a pass rate of 90% in three months of testing (20.04.2022 – 14.07.2022; 10 batches). Columbia Blood Agar Base (CM0331Y) has demonstrated a pass rate of 90% in eight months of testing (15.11.2021 – 04.07.2022; 10 batches). This shows that the performance is reproducible.

No data are available for the Columbia Blood Agar Base (CM0331M) package variant at the time of file writing, since it constitutes a new variant intended to be marketed during 2023. As for all the other packaging variants of Columbia Blood Agar Base (CM0331B, CM0331K, CM0331R, CM0331T, CM0331U, CM0331W and CM0331Y) containers, Thermo Fisher Scientific will constantly monitor any QN data, and any product packaging-complaints related to Columbia Blood Agar Base (CM0331M). Any anomaly recorded for the Columbia Blood Agar Base (CM0331M) variant will be treated with utmost urgency and thoroughly investigated.

For target organisms, when using 10-100 cfu inoculum of *Staphylococcus aureus* ATCC® 25923™ and incubating the device 37°C for 24 hours, the user can recover organisms with colony size and morphology as listed in this document. For target organisms, when using 10-100 cfu inoculum of *Clostridium sporogenes* ATCC® 19404™, *Clostridium sporogenes* ATCC® 11437™ and incubating the device 37°C for 48 hours under anaerobic conditions, the user can recover organisms with colony size and morphology as listed in this document. For target organisms, when using 10-100 cfu inoculum of *Streptococcus pyogenes* ATCC® 19615™, *Streptococcus pneumoniae* ATCC® 6305™ in media enriched with 5% v/v horse blood after incubation at 37°C for 24 hours, the user can recover organisms with colony size and morphology as listed in this document.

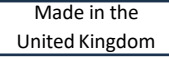
For target organisms, when using 10-100 cfu inoculum of *Neisseria gonorrhoea* NCTC 11148 in media enriched with 5% v/v horse blood and after incubating the device after incubating the Columbia Blood Agar Base (CM331) devices at 37°C for 48 hours under microaerophilic conditions, the user can recover organisms with colony size and morphology as listed in the document. For target organisms, when using the specified inoculum *Haemophilus influenzae* ATCC® 49247™ and *Haemophilus parainfluenzae* ATCC® 33392™ on plain media at 37°C for 18 hours, the user can measure zones of growth as listed in this document. For target organisms, when using the specified inoculum *Streptococcus pyogenes* ATCC® 19615™ on media enriched with 7% v/v horse blood plates and incubated at 37°C for 18 hours, the user can measure zones of growth as listed in this document. For target organisms, when using 50-120 cfu inoculum, of *Campylobacter jejuni* ATCC® 29428™, *Campylobacter jejuni* ATCC® 33291™ and *Campylobacter coli* ATCC® 43478™ and incubating the device at 41.5 ± 2°C for 44 ± 4 hours, the user can recover organisms with colony size and morphology as listed in this document. For target organisms, when using 100 to 10000 cfu inoculum, of *Streptococcus pyogenes* ATCC® 19615™, *Streptococcus pneumoniae* ATCC® 6305™, *Staphylococcus aureus* ATCC® 25923™ and *Escherichia coli* ATCC® 25922™ when tested in accordance with CLSI M22 A (incubation at 35°C for 18-24 hours), the user can recover organisms with colony size and morphology as listed in this document.

Bibliography

- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. "Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. "Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

Symbol Legend

Symbol	Definition
	Catalogue number
	In Vitro Diagnostic Medical Device
	Batch code
	Temperature limit

	Use-by date
	Keep away from sunlight
	Consult instructions for use or consult electronic instructions for use
	Contains sufficient for <n> tests
	Do not use if packaging damaged and consult instructions for use
	Manufacturer
	Authorized representative in the European Community/European Union
	European Conformity Assessment
	UK Conformity Assessment
	Unique device identifier
	Importer - To indicate the entity importing the medical device into the locale. Applicable to the European Union
	Made in the United Kingdom

©2022 Thermo Fisher Scientific Inc. All rights reserved. ATCC and ATCC catalogue marks are a trademark of American Type Culture Collection.
NCTC and NCTC catalogue marks are a trademark of National Collection of Type Cultures.
All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

For technical assistance please contact your local distributor.

Revision information

Version	Date of issue and modifications introduced
4.0	2023-11-10

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1. Public Health England. 2021a. „Identification of Streptococcus species, Enterococcus species and morphologically similar organisms.” UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. „Identification of Staphylococcus Species, Micrococcus Species and Rothia Species”. UK Standards for Microbiology Investigations.ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. „identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. „идентификация на Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Literatura

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations.ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. 'identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Litteratur

1. Public Health England. 2021a. "Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Tilgængeligt 15. oktober 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. "Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20". Vol. 14.
3. Public Health England 2020b. "Identification of Staphylococcus Species, Micrococcus Species and Rothia Species". UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. "Identification of Neisseria species. ID6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. "Identification of Clostridium species. ID8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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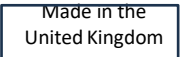
Bibliographie

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Zugriff am 15. September 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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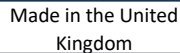
- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Τεύχος 4. Πρόσβαση 15 Οκτωβρίου 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Oxoid Limited, Wade Road, Basingstoke,
RG24 8PW, H.B.

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1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Último acceso: 15 de octubre de 2021. <https://www.gov.uk/government/publications/smi-id-4-identification-of-streptococcus-species-enterococcus-species-and-morphologically-similar-organisms>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>
4. Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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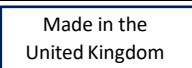
Bibliographie

- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. "Identification of Neisseria species. ID6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. "Identification of Clostridium species. ID8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Oxoid Limited, Wade Road, Basingstoke,
RG24 8PW, Royaume-Uni

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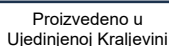
Bibliografija

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Pristupljeno: 15. listopada 2021. <https://www.gov.uk/guidance/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>
4. Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Hozzáférés: 2021. október 15. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>
- Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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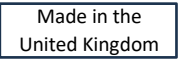
Bibliografia

- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Consultato il 15 October 2021. <https://www.gov.uk/government/publications/smi-id-4-identification-of-streptococcus-species-enterococcus-species-and-morphologically-similar-organisms>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>
- Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Oxoid Limited, Wade Road, Basingstoke,
RG24 8PW, UK

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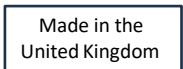
Literatūra

- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Prieigos data: 2021 m. spalio 15 d. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. 'Identification of Neisseria species. ID6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. 'Identification of Clostridium species. ID8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Oxoid Limited, Wade Road, Basingstoke,
RG24 8PW, JK

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Bibliografi

- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/government/publications/smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Bibliografia

- Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
- Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
- Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations.ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
- Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
- Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Bibliografia

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Edição 4. Acesso a 15 de outubro de 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. 'identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Bibliografie

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accesat pe 15 octombrie 2021. <https://www.gov.uk/government/publications/smi-id-4-identification-of-streptococcus-species-enterococcus-species-and-morphologically-similar-organisms>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations.ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. "Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. "Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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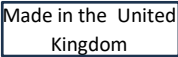
Zdroje

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. 'Identification of Neisseria species. ID 6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'Identification of Clostridium species. ID 8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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- Campylobacter Selective Supplement (Skirrow) (SR0069E/H)
- Staph/Strep Selective Supplement (SR0070E)
- Campylobacter Selective Supplement (Butzler) (SR0085E)
- Campylobacter Selective Supplement (Blaser-Wang) (SR0098E)
- Streptococcus Selective Supplement (SR0126E)
- Helicobacter pylori Selective Supplement (DENT)(SR0147E)
- Modified CNA Selective Supplement (SR0176G)
- Campylobacter Growth Supplement (SR0232E)

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Bibliografi

1. Public Health England. 2021a. 'Identification of Streptococcus species, Enterococcus species and morphologically similar organisms. UK Standards for Microbiology Investigations. ID 4 Issue 4. Accessed 15 October 2021. <https://www.gov.uk/uk-standards-for-microbiology-investigations-smi-quality-and-consistency-in-clinical-laboratories>
2. Public Health England 2020a. 'Group A Streptococcal Infections: Third Report on Seasonal Activity, 2019/20'. Vol. 14.
3. Public Health England 2020b. 'Identification of Staphylococcus Species, Micrococcus Species and Rothia Species'. UK Standards for Microbiology Investigations. ID 7 (4). <https://www.gov.uk/government/publications/smi-id-7-identification-of-staphylococcus-species-micrococcus-species-and-rothia-species>.
4. Public Health England. 2015. 'Identification of Neisseria species. ID6 <https://www.gov.uk/government/publications/smi-id-6-identification-of-neisseria-species>
5. Public Health England. 2016. 'Identification of Clostridium species. ID8. <https://www.gov.uk/government/publications/smi-id-8-identification-of-clostridium-species>

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Alla övriga varumärken tillhör Thermo Fisher Scientific Inc. och dess dotterbolag.