

CD49d (Integrin alpha 4) Monoclonal Antibody (9F10), PE, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human, Rhesus monkey
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PE, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	9F10
Conjugate	PE
Excitation/Emission Max	565/576 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_10717245

Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	5 µL (0.125 µg)/test	16 Publications

Product Specific Information


Description: The 9F10 monoclonal antibody reacts with human CD49d, the 150 kDa integrin alpha4 subunit. The complex of CD49d non-covalently associated with integrin beta1 (CD29), also known as VLA-4, is a receptor for fibronectin and VCAM-1 (CD106). This complex is expressed by thymocytes, peripheral lymphocytes, monocytes and eosinophils. CD49d also associates with integrin beta7 and binds to the Mucosal Addressin Cell-Adhesion Molecule-1 (MadCAM-1).^M

Applications Reported: The 9F10 antibody has been reported for use in flow cytometric analysis.^M

Applications Tested: This 9F10 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.125 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.^M

Excitation: 488-561 nm; Emission: 578 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.^M

Filtration: 0.2 µm post-manufacturing filtered.

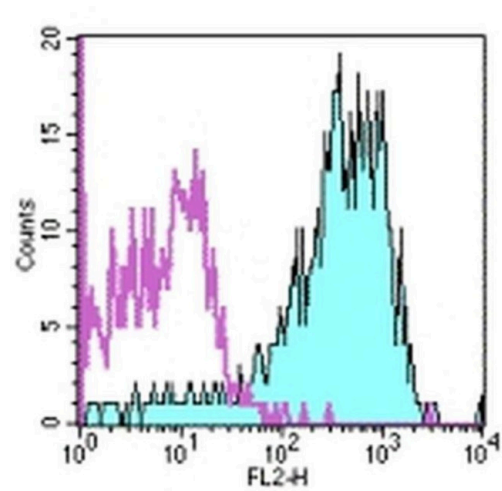

CD49d (Integrin alpha 4) Antibody (12-0499-42) in Flow
Staining of normal human peripheral blood cells with Mouse IgG1 kappa Isotype Control PE (Product # 12-4714-81) (open histogram) or Anti-Human CD49d (Integrin alpha 4) PE (filled histogram). Cells in the lymphocyte gate were used for analysis.

17 References

Immunocytochemistry (1)

<p>PloS one</p> <p>Tetraspanins CD81 and CD82 facilitate 41-mediated adhesion of human erythroblasts to vascular cell adhesion molecule-1.</p> <p>"12-0499 was used in Flow cytometry/Cell sorting to investigate whether any tetraspanin could modulate integrin function, showing that tetraspanins CD81 and CD82 facilitate 41-mediated adhesion of human erythroblasts to VCAM-1."</p> <p>Authors: Spring FA,Griffiths RE,Mankelow TJ,Agnew C,Parsons SF,Chasis JA,Anstee DJ</p>	<p>Year 2013</p> <p>Species Human</p>
---	---

Flow Cytometry (16)

<p>Stem cell research & therapy</p> <p>Mesenchymal stem cell aggregation mediated by integrin 4/VCAM-1 after intrathecal transplantation in MCAO rats.</p> <p>"Published figure using CD49d (Integrin alpha 4) monoclonal antibody (Product # 12-0499-42) in Flow Cytometry"</p> <p>Authors: Ran Y,Dong Y,Li Y,Xie J,Zeng S,Liang C,Dai W,Tang W,Wu Y,Yu S</p>	<p>Year 2022</p>
---	----------------------

<p>Communications biology</p> <p>Identification of potential chemical compounds enhancing generation of enucleated cells from immortalized human erythroid cell lines.</p> <p>"Published figure using CD49d (Integrin alpha 4) monoclonal antibody (Product # 12-0499-42) in Flow Cytometry"</p> <p>Authors: Soboleva S,Kurita R,Ek F,Åkerstrand H,Silvério-Alves R,Olsson R,Nakamura Y,Miharada K</p>	<p>Year 2021</p>
---	----------------------

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (i) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (ii) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (iii) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (iv) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.